

beauty

wellness

podo



BENTLON

English

BENTLON



Bentlon Beauty

Treatment chairs	8
Stools	18
Work lamps	20
Instrument tables	22
Softsteamers	24
Bentlon IPL	28
Bentlon DermaPerfect	32
Bentlon Micrabrader	36
Beauty equipment	40



Bentlon Wellness

Wellness couches	48
Wellness equipment	52



Bentlon Podo

Podiatry chairs	58
Podiatry devices	66
Podiatry S Cabinet	70
Cryopen by Bentlon	72
Podiatry accessories	74



Ramon Spitters
President Dancohr Corporation BV
Business Development Manager

Bentlon stays innovative!

With pride we present you the new Bentlon brochure which shows all new innovations and high quality products within the Bentlon assortment.

Bentlon products are developed and produced by Dancohr Corporation BV in Weert, the Netherlands. With 40 years of experience in the branch and an own R&D department, Bentlon is a true statement in the beauty, podiatry and wellness market. Mr. Ramon Spitters is the second generation of the Spitters family that leads the general management since 2012. He is also the driving force behind all new Bentlon innovations. As a family business by origin, personal contact with customers in the Netherlands and abroad is greatly valued.

For detailed product information, you can also browse the website on www.bentlon.com or find a local Bentlon distributor for personal advice. The possibility exists that not all Bentlon products in this brochure are available in your country. Please contact your local distributor for more information.

To contact us directly you can mail to info@dancohr.com or call +31 495 434 612. We hope you find great pleasure and inspiration going through this new edition.

Kind regards,
Dancohr Corporation.



Touch of Excellence

Well-balanced product collection

Bentlon treatment chairs are available in many different models with different features to ensure the right chair is available for anyone. From entry-level models with partly electric control to very advanced models with all kinds of characteristics that make your work easier. In terms of changing position of the chair, Bentlon offers great ease of use. Height, back, leg, and tilting position can be adjusted with all chairs, either electrical with 2, 3 or 4 actuators or with a pull-over ride. Standard models of the treatment chairs are often extendable with extra options such as a foot control, wheels, other head rests and 37 colors upholstery. The Bentlon assortment offers a solution for everyone!

Tuned to the market

A lot of research is executed before developing a new Bentlon chair. All new developments are carefully tuned to market trends and market needs. This is a must for Bentlon. Interesting suggestions from the market are always taken in serious consideration when starting new projects. Our R&D department works closely together with experts in different fields of practice such as our customer beauticians, podiatrists, physical therapists and others.

Excellent quality by experience

Functionality and user friendliness combined with a trendy design are leading for Bentlon furnishing and equipment. These aspects are the main criteria in a development project in which our customers take a central place. Bentlon strives to introduce new standards in the market which enhance your work experience, increase the ease of use and therefore improve your working position. Results are the patented 3-dimensional head rest, massage head rest, step-out assist and Lower Arm Support.

Quality and warranty

The high quality of Bentlon can be guaranteed by the fully computerized production processes with robots and constant quality control. Each product undergoes at least 20 quality checks during the production process which is also visible in the detailed and secured finish. Because of this, Dancohr gives 2 years general warranty and even 5 years on production errors.



10 reasons for a Bentlon choice

1. Quality

With 40 years of experiences in the branch Bentlon is the benchmark for quality.

2. Complete assortment

Bentlon offers a complete assortment for beauty, wellness and podiatry

3. Dutch production

The production process takes place at the local factory facilities in the Netherlands. It is fully computerized with constant quality checks to ensure Bentlon quality.

4. ISO certified

All internal labour processes are executed and audited according to European standards.



5. For specialists, by specialists

During each development project input and know-how is gathered from industry experts. Whether it is a customer beautician or a podiatrist of art designer, Bentlon works closely together to reach the best outcome and develop products that meet our clients needs.

6. Innovative

Bentlon is market leader with high quality, trendy products. A lot of time and effort is invested in the development of annual product introductions.

7. Patents

Several Bentlon Innovations are patented, such as the Lower Arm support, the unique Step Out Assist and the 3 Dimensional Head read for the beauty treatment chairs.

8. Quick Return on investment

Quick return on investment because of good price quality ratio. After gaining back the investment of a Bentlon Specials beauty device, the device will stay a profitable device for your salon!

9. Worldwide distribution

Bentlon is constantly focused on growing brand awareness worldwide. Currently the Bentlon products are exported to over 60 countries across the globe..

10. Excellent support

Bentlon offers a total concept including extensive coaching, training and support.



BENTLON

Beauty



Touch of Excellence



Treatment chairs	8
Stools	18
Work lamps	20
Instrument tables	22
Softsteamers	24
Bentlon IPL	28
Bentlon DermaPerfect	32
Bentlon Micrabrader	36
Beauty equipment	40



BENTLON



Beauty



Ref. 60006.100

...without compromising Bentlon quality!

- Functional design at entry-level pricing
- Optimal ease of use because of 2 actuators
- High working position
- Lifting capacity of 200 kg
- Innovative one touch treatment positions
- Ergonomic high quality upholstery





The Bentlon Bronze treatment chairs are equipped with 2 actuators which make it possible to electrically adjust the height and tilting function of the chair to set it to the desired position. It is also possible to put the chair in massage position electrically. The back rest adjustment is easy to operate by means of a gas compression spring. This combination of adjustment possibilities offers the specialist exceptional ease of use.



The treatment position of the Bentlon Bronze ensures a comfortable seated position for the client at all times. Thanks to the narrow back rest, the chair has a slender appearance, which makes the salon or consultation space appear more spacious. In addition, the arm supports can fold up completely, which makes it easy for the client to step into and out of the chair from the side.

To give the Bentlon Bronze chairs a luxurious appearance, Bentlon has selected high-quality imitation leather with a calfskin grain and a subtle matte-lustre finish. The upholstery is very easy to clean and is resistant to disinfectants. The abrasion- and tear-resistant upholstery ensures a longer service life.



Clients can step into and out of the chair easily by using the unique, patented electric step-in and out seat assist. This step-in and out seat assist lowers the chair, tilts the seat forward and rotates the leg and foot support inward. The chair follows and supports the client's movement as they step in or out. The chair also features robust construction which significantly increases the lifting and tilting power.



Ref. 60006100

Extra-low step-out seat assist

The extra-low step-out seat assist makes it even easier to step into and out of the chair. The seat of the chair can be lowered to a height of 51 cm. For elderly or people who have trouble with mobility, this provides good support—only a little effort is required to step into and out of the chair.

High working position

The Bentlon Beauty Bronze is characterized by an exceptionally high working position. The seat of the chair can in fact be raised to a height of 85 cm, which also allows the specialist to give a massage at this working height. If the specialist wants to perform a foot treatment, then that is also possible. In that case, the working height can be as high as a maximum of 93 cm. This allows the specialist to comfortably work on the client's feet while standing.

Foot switch

It is also possible to operate the Bentlon Beauty Bronze using a foot switch. This makes it easy to change position during the treatment, while maintaining hygienic working conditions.

BENTLON



Beauty



Ref. 60045.100

The logical next step!

- Combines ease of use with excellent lifting capacity and comfort
- Electrical seat, back and tilting position because of 3 actuators
- Foot control with pre-defined treatment position
- Foldable arm rests to easily get in and out of the chair
- Modern design and available in 37 colors
- Ergonomically shaped high quality upholstery



By means of a foot control the Bentlon Beauty Silver can be set to almost any desired treatment position. The foot control also has a home-position to easily bring the chair back to the base position. This means that almost all positions can be set using the electric foot control. The leg support can be tuned to the right angle using a pull-override. Also for clients getting in and out of the chair is easy because of the fold-away arm and leg support enabling them to get in and out of the chair from the side.



The Bentlon Beauty Silver has a very modern look and provides a very comfortable seating position to clients because of the ergonomically shaped upholstery which follows the contours of the body. The upholstery is made of high-quality leather and is available in 37 colors. The head support can be pulled out and can be adjusted to set the most ideal treatment position for any client.

Accessoires and optionals

The Bentlon Beauty Silver can be extended with optionals such as a seat cover, paper roll holder, wheels and a manicure tray. Moreover, you can choose to mount a massage head support and a Lower Arm Support.



Massage head support

When desired a C-curved head support, ideal for massages, can be mounted to the Bentlon Beauty Silver. It supports the face during a treatment in supine position. It offers the client extra comfort because of the breathing hole enabling him or her to optimally enjoy the treatment.

Lower Arm Support

To support your arms during a foot treatment it is possible to mount a Lower Arm Support which unburdens your shoulders and backs and gives you a better working position.

BENTLON



Beauty



Ref. 60025.100

Proven comfort and quality

- Stable high quality chair with 3 actuators
- Unique patented Step Out - assist
- Lifting capacity of 160 kg.
- Modern design and optimal seating comfort
- Imitation leather upholstery available in 37 colors!
- Also available with thermo/vibra functionality



The Bentlon Gold has the possibility to electrically set the right position of the height, back and tilting angle. The leg support can be tuned to the right position using the pull-override. Besides that, clients can get in and out of the chair effortless because of the fold-away arm and leg support which makes this chair very practical in use.



The Bentlon Gold has a very modern look and provides a very comfortable seating position to clients because of the ergonomically shaped upholstery which follows the contours of the body. The upholstery is made of high-quality leather and is available in 37 colors. This variation in combination with a colored foot frame (white, wengé or aluminium) makes it possible to customize the chair exactly to your wishes!



Clients can step into and out of the chair easily by using the unique, patented electric step-in and out assist. This step-in and out assist lowers the chair, tilts the seat forward and rotates the leg and foot support inward. The chair follows and supports the client's movement as they step in or out. The chair also features robust construction which significantly increases the lifting and tilting power.

Bentlon chairs can be fit-ted with relax vibration and chair heating. The chair heating is not only especially pleasant for the client, but is also energy saving because the temperature in the salon can be set a few degrees lower.



3-Dimensional adjustable head support

Also the Bentlon Gold is equipped with a 3-dimensional adjustable head support which makes it possible to not only move the head support higher but also to the back or frontwards, resulting in more comfort during a treatment in flat position!

Accessoires and optionals

The Bentlon Beauty Silver can be extended with optionals such as a seat cover, paper roll holder, wheels and a manicure tray. Moreover, you can choose to mount a massage head support and a Lower Arm Support.



BENTLON



Beauty



Ref. 60050.100

Beauty bestseller!

- Full electric control because of 4 actuators
- Foot control with home position, 2 memory positions and automatic massage position
- Solid chair with lifting capacity of 225 kg
- Luxurious edition with heating and massage
 - Heating system with 3 zones
 - New massage system with 5 modes and 2 programs



The Bentlon Platinum makes a statement in the beauty market. This stable and high-quality chair offers unique ease of use. Full electric control by means of a foot control enables the professional beautician to easily set the chair in the right position with the least of effort. With just one push on the foot control all actuators can be activated to put the chair in the desired treatment position, to activate the patented step out-assist and to tune the position of the leg support.



Bentlon chairs can be fitted with relax vibration and chair heating. The chair heating is not only especially pleasant for the client, but is also energy saving because the temperature in the salon can be set a few degrees lower.

This also creates for you as beauty specialist a considerably more comfortable working climate. Many other expansion options are available for your treatment chair, such as a paper roll holder, castors, a mountable manicure tray or the Lower Arm Support.



The new chair has a stylish and trendy design and its ergonomically shaped upholstery ensures that your client sits comfortably. This upholstery is made from high-grade imitation leather and is resistant to cleaning products without any solvents.



The unique electric step-in and step-out assist helps your clients get into and out of the chair without any trouble. This step-in and step-out assist lowers the chair, tilts the seat forward and turns the leg and foot section inwards. This follows and supports the client's stepping in and out movement. Moreover, the chair has a solid build, which considerably increases the ability to hoist and lift.

3-Dimensional adjustable head support

Naturally, this multifunctional patented head support can be raised, brought down and tilted, but it also offers the option to adjust it backward or forward. This offers advantages when the client lies on the back or the stomach, because your client's neck is stretched better for a treatment and your client is in a more relaxed position.



Massage head support

With certain treatments the use of the special C-curved massage head support is preferred. For treatments in prone position the head support gives complete support around the face, resulting in greater comfort during a longer treatment. The length and depth of the head support can also be adjusted.

Lower Arm Support

The media both at home and abroad has been full of the Lower Arm Support, a special support for pedicures that supports the lower arms during the treatment so that the load is taken off the shoulders and back and an improved working posture is assumed almost automatically.

Lower Arm Support



Together with the Delft University of Technology, Bentlon has made an important invention with a worldwide patent: the Lower Arm Support. This is an ergonomic support with bin for podiatrists and represents a breakthrough in the podiatry industry. The work of the podiatrist requires a high level of precision, concentration and long periods of muscle contraction. During the treatment a certain posture is adopted or the load on the body is on one side. Often podiatrists sit with a bent posture or in a slightly twisted position to obtain a good view when working and their arms are continuously raised while working. This is very tiring and wearing for the body.

Support

The Lower Arm Support supports the lower arm so that the weight of the upper arm, lower arm and hand does not have to be supported wholly by the muscles in the back and shoulders. Because of this support the treatment chair that the patient sits on can be set higher; this means that the podiatrist can see the foot better and is able to adopt a more upright position of the lower back and the neck. This means a considerable improvement in the work posture and is much less tiring for the shoulders and arms.

Hygienic

The Lower Arm Support also has a removable bin that is easy to remove and put back and is easy to clean due to its funnel shape. The support bracket and the bin have been developed so that there is enough space around the foot to work on it. The Lower Arm Support is suitable for left-handed and right-handed podiatrists. The Lower Arm Support can be fitted on all Bentlon podochairs. You can also attend a workshop that will teach you how to make optimum use of the Lower Arm Support so that you considerably improve your work posture and can save your body from painful complaints due to an incorrect position.



Treatment chairs	Bronze	Silver	Gold	Platinum
Ref.no. with hand control	60006	60045	⊗	⊗
Ref.no. with foot control	60008	⊗	60025	60050
Ref.no. with chair heating	⊗	⊗	60026	60051
Ref.no. with vibration and heating	⊗	⊗	⊗	60052
Number of actuators	2	3	3	4
Length	185 cm	178-211 cm	195-225 cm	178-211 cm
Chair width	75 cm	80 cm	79 cm	80 cm
Sitting width	55 cm	80 cm	61 cm	57 cm
Height	51-83 cm	60-94 cm	65-89 cm	60-94 cm
Chair weight	53 kg	60-65 kg	54-60 kg	60-65 kg
Max. permitted weight	200 kg	225 kg	160 kg	225 kg
Watt	96W	230W	230W	230W
Upholstery	Imitation leather	Imitation leather	Imitation leather	Imitation leather
Coloured upholstery	1	37	37	37
Aluminium	⊗	⊗	○	○
Wengé	⊗	⊗	○	○
Extendable head support	●	●	●	●
Tilting head support	⊗	⊗	●	●
Head support depth adjustment	⊗	⊗	●	●
Extendable foot support	⊗	●	●	●
Fold-away arm supports	●	●	●	●
Memory lying and exit position	⊗	⊗	●	●
Height adjustment	Electrical	Electrical	Electrical	Electrical
Tilt adjustment	Electrical	Electrical	Electrical	Electrical
Back adjustment	gasveer	Electrical	Electrical	Electrical
Leg adjustment	Gas compression spring	Gas compression spring	Gas compression spring	Electrical
Link leg and backrest	⊗	⊗	Mechanical	Electrical
Patented step-out assist	●	⊗	⊗	●

For all the colors of the upholstery go to page 76

Accessories	Ref.
Chair covering towelling	61040/41/42
Protective set towelling	61045
Manicure tray with brace	61080 + 61079
Paper roll holder with paper roll	61087
2-ply paper roll 50cm x 50cm	61084/61085
Lower Arm Support Beauty Blue	61088.100
Lower Arm Support Beauty White	61853.100
Support roll round	61090
Support roll semi-circular	61091
Wheels	61099
Massage head support	61083
Cover Massage head support	61022



Ref. 61088



Ref. 61080 + 61079



Ref. 61087 + 61085

BENTLON

Stools



Ref. 60160

For the correct working position

- Varied assortment
- Ergonomic upholstery
- In height adjustable
- With or without back rest
- In matching color with treatment chair!

Choosing the right stool

Stools are available in various designs. The stool is very important for correct work posture and should be comfortable. This is a purely personal thing. In the Bentlon stool range you will find exactly the right stool for you. In addition to a wide selection of cover colours, with most models you can choose from a neutral grey or chrome base. All stools are stable and have five legs with castors. The height adjustment is continuous with a hydraulic gas compression spring. On most models the back support can be adjusted both in height and depth so that you can always set the perfect sitting position.



Stools	Bronze	Silver Comfort	Silver Plus	Silver Pro	Silver Ergo	Gold	Platinum	Platinum+
Ref.no.	50145	50155	50165	50175	50190	60150	60155	60160
Seat	Round	Round	Saddle	Oval	Pony saddle	Round	Round	Oval
Seat width	38 cm	38 cm	38 cm	42 cm	59 cm	38 cm	38 cm	45 cm
Height	48 - 61 cm	48 - 61 cm	48 - 61 cm	48 - 61 cm	54 - 72 cm	43 - 55 cm	43 - 55 cm	46 - 61 cm
Back support	⊗	●	●	●	⊗	●	●	●
Adjustability back sup.	⊗	Height+depth	Height+depth	Height+depth	⊗	Depth	Depth	Height+depth
Max. permitted weight	120 kg	120 kg	120 kg	120 kg	120 kg	120 kg	120 kg	120 kg
Rotation	360°	360°	360°	360°	360°	360°	360°	360°
Base	Chrome	Chrome	Chrome	Chrome	Chrome	Chrome	Chrome	Chrome
Base width	59 cm	59 cm	59 cm	59 cm	59 cm	59 cm	59 cm	59 cm
Castors	5	5	5	5	5	5	5	5
Weight	5,5 kg	7,8 kg	7,8 kg	7,8 kg	8,5 kg	11,8 kg	11,8 kg	11,8 kg
Upholstery	Imitation leather	Imitation leather	Imitation leather	Imitation leather	Imitation leather	Imitation leather	Imitation leather	Imitation leather
Coloured upholstery	37 colours	37 colours	37 colours	37 colours	37 colours	37 colours	37 colours	37 colours
Aluminium	⊗	○	⊗	⊗	⊗	⊗	⊗	○
Wengé	⊗	○	⊗	⊗	⊗	⊗	⊗	○
Adjustable seating	⊗	⊗	⊗	⊗	⊗	⊗	○	⊗

For all the colors of the upholstery go to page 76

Accessories	Ref.
Stool covering toweling Silver/Gold	51150
Bin tabouret Platinum+	61105
Bottle holder tabouret Platinum+	61106

Ref. 61105



Ref. 61106



BENTLON

Work lamps



Ref. 50703

To have the perfect visibility!

- Professional work lamps
- Enlarge to minimal 3x
- Varied assortment
- High quality because of glass, lens, light and construction

A work lamp remains to the basic equipment of any beauty salon. Its quality is determined by the quality of the glass, the lens, the light and the construction of the lamp. The most advanced work lamps is the best choice for longterm use which needs good light, good color display of the skin and no distorted view. If a work lamp is chosen purely for extra light, then a basic version will do.

Bentlon Bronze

This light work lamp is easy to use and does not emit heat. With 141 cm. this lamp has the longest working arm of the range which makes it possible to work from a distance. The Bentlon Bronze work lamp also has a so-called Free Motion system which enables you to easily place the lamp to the desired position. The lamp enlarges up to 3x and contains a 28w high frequent TL light which gives a very clear undistorted view.

Bentlon Silver

What makes this work lamp special is the double lighting from the side which results in a shadow-free light. The lights can also be used separately from each other. The glass enlarges up to 3x.

Bentlon Wave Plus

This rectangular, broad spectrum work lamp gives a very comfortable and clear view. The Wave Plus also has double lighting from the side for a shadow-free light.

Bentlon Circus

This work lamp has an operating arm of 100 cm and a round lens which is provided with a protection cap. The T5 22w TL lighting offers a clear and undistorted view.

Work lamps	Bronze	Silver	Wave Plus	Circus
Ref. no.	50703	50708	50721	50725
Length arm	141 cm	132 cm	105	100 cm
Lens dimension	Ø 17.5 cm	19x17.5 cm	22,7x22,4	Ø 16.5 cm
Lens magnification	3x	3x	3,5x	3,5x
Light power	28 W	2x 13 W	2x 6 W	22 W
Power supply	230 V	230 V	18 V	230 V
Table clamp	●	●	●	●

Accessories	Ref.
Ring lamp 22W	51702
Table clamp Wave II Plus	51719



Ref. 50703



Ref. 50708



Ref. 50721



Ref. 50725

Instrument tables and cabinets



Ref. 65625.06

Ref. 65625.00

Practical interior solutions

- Practical in use
- Extendable with UV drawer or tableau
- Optional work lamp holder
- Various designs

Instrument tables

An instrument table is vital in the salon. Bentlon instrument tables are constructed in such a way that you can use them whilst having the products you need for treatment close at hand, but the instrument tables are also suitable for placing equipment in or on. With each model you have the choice of two or three shelves. There is the option to have a drawer under each shelf to store small materials. All instrument tables have damped castors. You can place a magnifying lamp on the table top with a special holder.



Ref. 60300

Ref. 60305

Bentlon Gold Cabinet

Have you recently purchased a Bentlon Wellness Couch Gold? And would you like to give your salon an even more coordinated look? If so, the Bentlon Gold Cabinet may be just what you are looking for, since it is available in the same colours as the Bentlon Wellness Couch Gold. This cabinet offers a highly stable and mobile work surface for placing your equipment and accessories. It has three spacious drawers and the top, which is made of hardened glass and therefore hygienic and easy to clean, offers ample space for accessories and products.



Ref. 65625.00

Ref. 65625.01

Ref. 65625.02

Ref. 65625.03

Ref. 65625.04

Ref. 65625.05

Ref. 65625.06

Cabinets	Silver 2	Gold 3	Gold	Gold-uv-lade	Gold-tableau	Gold-tableau /uv lade
Ref.no.	60300	60305	65625	65626	65627	65628
Width	65 cm	65 cm	60 cm	60 cm	60 cm	60 cm
Depth	42 cm	42 cm	59 cm	59 cm	59 cm	59 cm
Height	69 cm	91 cm	80 cm	80 cm	80 cm	80 cm
Levels	2	3	1	1	1	1
Space between levels	40 cm	27,5 cm	⊗	⊗	⊗	⊗
Drawers	○ max. 2	○ max. 3	3	3	3	3
Wheels	●	●	●	●	●	●
Holder for magnifying lamp	○ 61305	○ 61305	●	○	●	○
Tray dimensions	⊗	⊗	⊗	⊗	34,5 x 36,5 cm	⊗

Accessories	Ref.
Drawer for Silver/Gold	61300
Holder for magnifying lamp	61305
Lid softsteamer tube for Silver/Gold	61310

Softsteamers



Modern design, practical and complete!

- Modern design
- Steam function adjustable from 1 up to 9 steam strengths depending on model
- Optimal ease of use because of rotatable arm
- Practical in use as a side table
- Work perfectly with tap water
- Optional socket, storage and work lamp holder

Bentlon introduces four new models of SoftSteamers all of which fulfil the needs of every aesthetic specialist. From the starter's model with its basic steam functionality up to the advanced steamer with a range of extras.

Complete range

Together, the new Bentlon SoftSteamer models make up a complete range in which each model offers additional functionality or options. This provides every beautician with the possibility to select a Bentlon SoftSteamer with the options that best fit her salon.

Softsteamer Bronze

Starter's model of superior quality! The Bentlon SoftSteamer Bronze is the starter's model steam device with authentic steam function. This model is equipped with a felt ring in which essential oils can be placed, which makes it possible to add an extra dimension to the steaming with the addition of aromas. The Bentlon SoftSteamer Bronze is equipped with a rotating steam arm. This arm makes it very easy to set the steam head and steam direction in the correct position without the need to adjust the height of the SoftSteamer itself. The modern design of the new Bentlon SoftSteamer Bronze, combined with its practical attachable top tray, make this SoftSteamer the ideal starter's model.

SoftSteamer Silver

Basic model with the possibility for variety! The Bentlon SoftSteamer Silver is a basic model of steam device which in addition to the basic steam function features a felt ring which gives the specialist the option for variety by adding herbs or using essential oils. In addition, this model also has an ozone function. The rotating steam arm on the Bentlon SoftSteamer Silver is not only very user-friendly and practical; it also offers the possibility of switching between two levels of steam strength. The design of the new Bentlon SoftSteamer The modern design of the new Bentlon SoftSteamer Silver, combined with the practical attachable top tray and the variable steam strengths, makes the Bentlon SoftSteamer Silver a solid choice for every salon.





Ref. 60290

SoftSteamer Gold

Top model with a range of options! The Bentlon SoftSteamer Gold is a top model steam device which offers a range of extra options in addition to the regular steam functionality. From a practical perspective, the SoftSteamer Gold is an ideal model—handsome finish, practical top tray, completely concealed power point and the option to add a universal magnifying lamp holder. This complete concept makes the Bentlon SoftSteamer Gold an excellent device for even the smallest spaces and makes it the pre-eminent must-have addition to every salon. In addition, the Bentlon SoftSteamer Gold has an integrated power point to which other Bentlon Beauty devices can be connected and a universal magnifying lamp holder.

SoftSteamer Platinum

All-in-one model of superior quality! The Bentlon SoftSteamer Platinum is the queen of the SoftSteamers. A top model which combines steaming with additional functionality and options for the ultimate all-in-one model! The design of the new Bentlon SoftSteamer Platinum is sleek and modern but above all innovative! From a practical point of view, the Bentlon SoftSteamer Platinum is a true workhorse with its integrated storage drawers on the front, a practical tray on top, an integrated power point and a universal magnifying lamp holder. It is this complete concept that makes the Bentlon SoftSteamer Platinum the rightful queen of the SoftSteamers!

The design of the new Bentlon SoftSteamer Platinum is sleek and modern but above all innovative! From a practical point of view, the Bentlon SoftSteamer Platinum is a true workhorse with its integrated storage drawers on the front, a practical tray on top, an integrated power point and a universal magnifying lamp holder. It is this complete concept that makes the Bentlon SoftSteamer Platinum the rightful queen of the SoftSteamers!

Softsteamers	Bronze	Silver	Gold	Platinum
Ref. no.	60255	60265	60280	60290
Rotatable steamarm	●	●	●	●
Integrated tripod	●	●	●	●
Power	1200 W	1200 W	1200 W	1200 W
Reservoir capacity	1,75 ltr	1,75 ltr	1,75 ltr	1,75 ltr
Shockproof reservoir	●	●	●	●
Warm-up time	8	8	8	8
Steaming time	25	25	25	25
Steam strenghts	1	2	9	9
Covered steam arm	●	●	●	●
Aromas	●	●	●	●
Herbs	⊗	●	●	●
Ozon	⊗	●	●	●
Standby position	⊗	●	●	●
Boil-dry protection	●	●	●	●
LED water indication	⊗	⊗	●	●
Beauty apparatuur tableau	⊗	⊗	●	●
Integrated power point	⊗	⊗	●	●
Magnifying lamp holder	⊗	⊗	○	○
Storage	⊗	⊗	⊗	●

Accessoires	Ref.
Bentlon Herbs dry skin	45050
Bentlon Herbs sensitive skin	45051
Bentlon Herbs oily skin	45053
Bentlon Herbs set	45055
Essential oil felt ring	51290
Steam outlet + 2 felt rings	60210



Ref. 45050

BENTLON



IPL



Ref. 60655

Permanent hair removal and skin rejuvenation

- **Can be used for multiple treatments**
 - Permanent hair removal
 - Skin rejuvenation
 - Rosacea
 - Pigment spots
 - Acne
- **Life time of at least 60.000 flashes**



The Bentlon IPL Dual Power Flash Gold

IPL is the abbreviation for Intense Pulse Light. It is a strong light beam with a specific wavelength. The amount of energy from this light beam, in this case no less than 18 joule/cm², is absorbed by the hair or skin melanin (pigment). The Bentlon IPL has specific effects and applications depending on the wavelength. The Bentlon IPL Dual Power Flash Gold has a water-cooled pulse light with a very long lifespan, up to at least 60,000 pulses. This guaranteed lifespan ensures the force and thus the energy level remains consistent. Because of the water cooling, the pulse light does not get hot and this means that you do not have to stop during the treatment to allow the lamp to cool down. The operation of the Bentlon IPL Dual Power Flash Gold is also unique. No complicated settings of joule/cm², you just enter the hair colour, the hair thickness and the skin type in the menu and the Bentlon IPL sets the correct light power and pulse duration itself.

The treatments

The new Bentlon Dual Power Flash Gold can be used for a lot more than just hair removal as the unit has a considerably broader area of application. The numerous treatment options allow you to reach a larger target group, which eventually results in a high return and short payback period.

Permanent hair removal

Intense Pulse Light is a method for the definitive treatment of excess hair, or more precisely, unwanted hair. It is a strong light beam with a specific wavelength. The hairs are exposed to the pulse light for only a fraction of a second. An amount of light (=energy) is absorbed by the pigment in the hair.

The pulse light adapts itself to the thickness of the hair. The pigment heats up, transmits the heat to the root and burns it. If the temperature in the hair root reaches 70 degrees, the hair follicles that supply the hair root are seared shut so that the hair cannot grow back there. The handle is continuously cooled. Within a fraction of a second dozens of hair follicles can be treated at the same time. This makes it a very quick and effective method of hair removal for almost all skin types.



Skin rejuvenation

Over the years, the characteristics of ageing become more pronounced. Because of the reduction in elastin fibres, the skin loses its elasticity. As it gets older, the skin develops wrinkles, slackens and becomes drier. A treatment with Bentlon IPL Dual Power Flash Gold is for people with slack skin, atrophic skin, and skin with wrinkles and fine lines. Specific zones and wrinkles can also be treated rather than the whole face. With very specific IPL treatment, wrinkle depth is reduced and small lines even disappear completely.

Treatment of rosacea

Rosacea is caused by the fact that small blood vessels shine through the skin. This makes the skin red and the redness becomes even worse with temperature changes, sunbathing or the use of alcohol. The energy from the IPL lamp closes these blood vessels, causing the redness in the skin to disappear. The energy from the IPL treatment can close these red veins without affecting the skin. The light spectrum is such that only the red veins absorb the energy and are sealed off.



Treatment of pigmented areas

Pigmented areas and hyperpigmentation can be very disfiguring, in particular if they appear on the face. This very specific treatment focuses on the darker pigmented areas and makes them several tints lighter or removes them completely. For the treatment of pigment, the entire colour spectrum is needed. Dark pigments absorb all wavelengths of visible light. During the treatments, the pigment which is deeper in the skin is brought to the surface. After the first two or three treatments, the pigmented area will therefore become darker.

Acne treatment

Acne is successfully treated with the Bentlon IPL Dual Power Flash Gold. The treatment fights infection, which results in a considerable improvement of the skin. The active acne that is present disappears and the formation of new infections is minimised. This treatment is suitable for the face, back and chest. A timely treatment of the acne with the Bentlon IPL prevents the formation of unattractive scars. The heat that is introduced into the skin in a controlled and dosed manner has a positive effect on the sebaceous material. The sebaceous material becomes more fluid and also the bonds between cells become more flexible, reducing the amount of infections. In addition, the blue and green light of the colour spectrum have an adverse effect on the P.acnes bacteria, causing it to self-destruct.

The result

With the Bentlon IPL Dual Power Flash Gold treatments no artificial substances are injected into the skin, but the energy from the IPL pulse light takes effect in the skin. The skin repairs and improves itself with long-lasting results. The process takes place evenly and starts immediately with the first treatment.

Experience it yourself

In this brochure we have tried to present as complete a picture as possible of the Bentlon IPL Dual Power Flash Gold. However, we believe that you will only be really convinced if you are able to work with the Bentlon IPL Dual Power Flash Gold yourself. We should like to offer you the opportunity to get to know this quick and effective way of hair removal during a de-mostration in your salon.

Total concept

If you opt for the Bentlon IPL Dual Power Flash Gold, you are automatically opting for optimal service and support. The Bentlon IPL Dual Power Flash Gold is not only a unit, but a total concept with marketing support, Bentlon contactgel and eyewear protection. You receive the training course and accompanying manual completely free of charge. You are taught the theory and practice of the treatments. Of course, you will receive a certificate as proof that you have completed the course.

Bentlon IPL Dual Power Flash Gold pulse lights

The Bentlon IPL Dual Power Flash Gold comes with two types of pulse lights: one for the body (Hair Removal) and one for the face (Skin Rejuvenation). The Hair Removal pulse light is included as standard in the Bentlon IPL Dual Power Flash Gold's total package, together with a detailed training course and marketing material. The wavelength of this light is very important because it determines the penetration depth into the tissue. A wavelength of 610 to 950 nm is most suitable for hair removal. The shorter the wavelength, the higher the energy content. Below is an overview of the two pulse lights to be used with the Bentlon IPL Dual Power Flash Gold.

Patented hand piece

The Bentlon IPL Dual Power Flash Gold has patented hands pieces with automatic energy compensation that assure a constant level of light energy, namely 60,000 flashes throughout its life. The glass part of the flashbulbs is more effective thanks to the new design. The flashbulbs penetrate deeper in the skin because of the beveled edges of the glass.

Device	IPL Dual Power Flash Gold device
Ref. no. 230 V	60655
Power	up to 18 J/cm ²
Width x depth x height	47,5 x 49,9 x 40,5 cm
Weight	28 kg
Watt	900 W
Carrying case	No
Worktable cabinet	No
Training and certificate	Yes
Hydraulic cooling system	Yes
Applicators	HR (basic) SR (optional) up to 1 flash every 3 sec.

	Hair Removal	Skin Rejuvenation
Applicator	HR (<i>basic</i>)	SR (<i>optional</i>)
Output energy	12 to 18 J/cm ²	6 to 18 J/cm ²
Size of treatment head	7.5 cm ²	3.75 cm ²
Wavelength	610 nm - 1100 nm	475 nm - 1100 nm
Skin type	1 - 5	1 - 4

BENTLON

DermaPerfect

Ref. 60435

Perfects your skin!

- Removes cosmetic skin imperfections successfully, effectively and quick
- Local and goal oriented, tissue friendly
- Surrounding skin tissue will not be effected or damaged
- Treats various skin imperfections

The Bentlon DermaPerfect will successfully, effectively and quickly remove cosmetic skin irregularities such as wrinkles, fibromas, and tiny red veins. The Bentlon DermaPerfect is driven by a monopolar RF signal (Radio Frequency Signal). This 4mHz Radio Frequency Signal causes vibrations that concentrate on the tip of a tiny, non-insulated nickel rod. These vibrations create heat development in the skin, concentrated on the skin irregularities. This way, the skin irregularities in or on the skin are coagulated or removed. The Bentlon DermaPerfect treats the skin extremely locally and focuses on the particular spot, whereby the surrounding skin is not affected. That is why this treatment is very kind to the tissue.

Control circuit

The Bentlon DermaPerfect is a compact device with a clear display showing the different functions. It has the capacity to treat small areas and give instantly visible results. Throughout the treatment, safety is of the utmost importance. The menu contains a number of treatments. Simply select the treatment concerned, after which the Bentlon DermaPerfect will automatically use the correct settings and percentages. This way, you can start the treatment at the touch of a button, without having to perform any unnecessary actions. This procedure presents a very user-friendly and safe working method.

The treatments

In general, the treatments that make use of the Bentlon DermaPerfect are painless and short. However, the chance exists that a second treatment is necessary. This depends on a number of factors, including the size, the time a particular skin irregularity has already been present, and the skin's reaction.



The Bentlon DermaPerfect will treat various types of skin irregularities:

- Spider Naevus
- Teleangiectasias (rosacea)
- Blood vesicles
- Individual hairs
- Fibromas (tiny warts)
- Pigmented patches
- Wrinkles
- Age spots

The result

The Bentlon DermaPerfect is a highly effective device with a convincing effect. The monopolar Radio Frequency Signal will ensure quick, painless results within just a few minutes. Often, the results are visible immediately after the treatment. In some cases, the treated area may develop a small scab. Make sure to inform the client not to touch the scab – this may negatively affect the healing of the wound.



Home care

In order to maintain the condition of the treated skin area, it would be wise to advise the client to use Bentlon Care products. This way, the consumer can make sure that his or her skin is well cared for and will heal quickly.

Total concept

Upon purchasing the Bentlon DermaPerfect including case, a total concept is provided. This means that all products, promotion and trainings needed to perform treatments using the Bentlon DermaPerfect are comprised in the concept.

Including for example, the following:

- The operating pedal
- The needle holder with connection cable
- Three sets of F Care systems and Ballet disposable needles (130 needles) of 0.08 mm and 0.15 mm.
- Convenient, made-to-measure storage case
- The training, to the value of EUR 99, is offered free of charge upon purchase of the Bentlon DermaPerfect. During this training, you will not only learn about needle techniques, but you will also gain knowledge that will help you make the right diagnosis.
- Marketing material, such as consumer brochures and posters.
- A Bentlon Care package for skin care during and after the treatment.


The following products are included in the total concept:

- 45024 Bentlon Cleansing Milk
- 45026 Bentlon Calming Lotion
- 45008 Bentlon Skin Rebuilding Serum
- 45012 Bentlon Skin Rebuilding Cream
- 45019 Bentlon Depigmentation Serum
- 45020 Bentlon SPF Cream 50+

Profit

The Bentlon DermaPerfect has a high return on investment. It will pay for itself quickly, because the skin irregularities are treated easily and quickly. Thanks to the Bentlon DermaPerfect, you are able to perform many treatments. This also means that your target group is very wide. Using the total concept of the Bentlon DermaPerfect, you can perform 130 treatments. On average, you will have earned back your investment after 20 hours of treatment using the Bentlon DermaPerfect. Naturally, it depends on the treatment, which costs are involved.

The costs vary from EUR 25 to EUR 135. It is important to keep in mind that your client will value a treatment with a permanent result more highly than a manual treatment. After all, the skin is treated within minutes and shows the desired result. Which is, of course, what every client would wish for after a visit to the beautician.

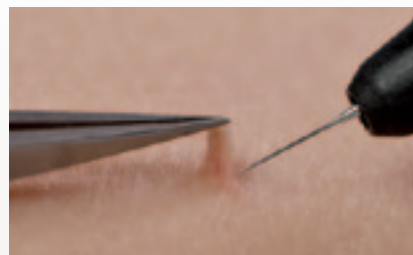

Ref. 61436



Ref. 61432 / 61433 / 61434



Ref. 61435

Device	DermaPerfect
Ref. no.	60435
System	Bentlon DermaPerfect
Output	Thermocoagulation signal
Width x depth x height	36 x 27,5 x 12 cm
Weight	10 kg
Power supply	Spanning 115-230V / 50-60 Hz
Power	115 VA
Wave frequency	4 mHz (+/- 0,1 mHz)
Max. output voltage	500 V
Max. duration of the impulse	0,8 s
FA/FB 500mA T (230V)	Reference melting fuse
FA/FB 1A T (115V)	Reference melting fuse
Display	Graphic display
Accessories	61432 K3i needles, 61433 K6i needles
	61434 R3i needles, 61435 case
	61436 Footswitch
	61437 needle holder with needle connection
	61438 needle holder

BENTLON



Micrabrader



Effective skin improvement for an even and revitalized complexion!

- Ideal treatment for skin rejuvenation and skin improvement
- Reduces skin imperfections and wrinkles to a minimum
- Results visible after the first treatment!
- High return on investment



Revitalize the skin and improve the skin tone!

The Bentlon Micrabrader Gold can be used to reduce wrinkles and various skin imperfections to a minimum. Microdermabrasion is used most often to fade out wrinkles and fine lines, rejuvenating the complexion. A microdermabrasion treatment is also perfectly suitable to make other skin lesions less visible, such as (acne)scars, stretch marks and pigmentation spots.

Microdermabrasion

By means of a special suction technique sterile microcrystals are blasted onto the skin, loosening old and dead skin cells which are drained with the microcrystals. The crystals remove the upper skin layer, stimulating the regeneration process of the skin at the same time. The regeneration capacity of the skin makes the skin heal itself by creating new and young cells. Besides that, blood circulation is stimulated during the procedure with the positive effect of activating the production of collagen and elastine fibers. Results are already visible after the first treatment. Results after an entire cycle are unique.

The Treatment

Microdermabrasion starts with cyclic treatments every 3 to 4 weeks until the desired result is reached. After that a treatment every 4 to 6 weeks is sufficient to keep good results. The treatment is very intensive, the client will feel the crystals flowing over the skin which is however not experienced as painful.



Bentlon Micrabrader Gold

The Bentlon Micrabrader Gold is the latest Bentlon development; a beauty device that offers excellent results but also is practical and hygienic in use and maintenance. Using soft touch keys, the beautician can choose the desired settings to start the treatment, either with soft or hard microcrystals. The Bentlon Micrabrader Gold has the special feature to memorize 4 preference settings to easily start a new treatment. Pressure of the crystal flow can also be adjusted and read from the digital display. A separate button will start the automatic cleansing function.

The crystal flow starts when the suction circuit is closed, simply when the hand piece is placed onto the skin. Once the hand piece is not in contact with the skin and thus the vacuum is interrupted, the crystal flow stops immediately. Because of this technique microcrystals will never be blown into the salon.



Bentlon Care products

Microdermabrasion is a very intensive treatment which requires good pre-and after care the client needs to abide by. The skin needs to be nourished with products that enable quick regeneration of the skin. Bentlon Care skin care products are precisely tailored to the aftercare of microdermabrasion.

During the treatment cycle the client should use Bentlon Care products at home to get the best results. After the cycle the client can return to using their regular skin care brand.



High return on investment

The Bentlon Micrabrader Gold has a very high return on investment. The treatment takes about as long as a regular facial, however the returns are remarkably higher. On average €85,-* can be charged per treatment of the entire face, entering a cycle of treating every 3-4 weeks.

** Average treatment prices can vary per country. Prices are indicative.*

Full support

When choosing the Bentlon Micrabrader Gold, intensive coaching and service is included. Moreover, various professional marketing materials are available to promote the treatment in your salon.

Warranty and maintenance

A Bentlon Micrabrader Gold is an investment with high returns for your salon. Bentlon recommends to maintain the device well and clean it on a regular base according to instructions in the manual. Good maintenance guarantees a longer lifetime.

All Bentlon equipment have a general warranty of 2 years and 5 years on construction errors.

Device	Micrabrader
Ref.no. 230 V	60491
Ref.no. 115 V	60493
Width x depth x height	37 x 37,5 x 23,5cm
Weight	10,5 kg
Watt	90 W
Treatment	Micro-dermabrasion
Operation	Soft-touch
Display	Digital
Power	< -0,6 bar
Memory positions	4
Automatic clean-function	●
Crystal	Aluminiumoxide (soft and hard)
Containers	2
Course and certificate	●

Accessories	Ref.
Microcrystals - hard	61502
Microcrystals - soft	61503
Starters pack (+/- 25 treatments)	61504
Pre-filter	61505
Input filter	61506
Nozzle set of 25	61507
Spray nozzle	61509



Ref. 61503



Ref. 61507

Beauty equipment



Ref. 60405



Ref. 60550



Ref. 60455

The perfect addition for your salon

- Brush device
- Combination suction/spray brush device
- Ionisation device
- High frequency device
- Cooling/heating device
- Compress heater
- Desinfection/sterilisation device

Brushing machine

The soft-touch keys on the completely flat display make the Bentlon brush machine very user-friendly. The brush can be set to rotate both clockwise and anti-clockwise and the speed can be adjusted accurately. The brushing machine is supplied with 1 pumice pad (4.5 cm) and 3 different brushes, namely a soft face brush (2 cm), a soft face/body brush (4.5 cm) and a hard body brush (4.5 cm). The set can be expanded with a cleaning sponge (4 cm) and an extra large body brush (6 cm). The holder for the pistol grip on the side and the two dishes on the unit mean that all loose accessories can be kept within reach.

Vacuum/spray/brushing machine

In addition to the features of the brushing machine described above, the Bentlon Combi machine also has a vacuum and spray function. The vacuum function is used to give anti-wrinkle treatments, to carry out pulsating massages and to release comedones in the skin. The spray function is equivalent to cooling vapour. To finish a treatment, an ultra-fine mist of a lotion can be used as a great way to refresh the skin.

Ionisation machine

Active ingredients can be introduced into the skin more effectively by means of ionisation, whereby a low galvanic current is used to pass in ampoules, aromas or other active ingredients. If the polarity of the ionisation machine is changed to desincrustation, blockages in the skin due to hardened sebaceous matter can be saponified, which allows deep cleansing to be performed. The ionisation machine is supplied with a roller/spherical electrode and a tweezer electrode. The digital display clearly shows the strength that has been set.

Device	Brush equipment	Combination unit	Ionisation unit
Ref.no. 230 V	60405	60455	60550
Width x depth x height	350 x 270 x 110 mm	350 x 270 x 110 mm	350 x 270 x 110 mm
Weight	4,2 kg	6,3 kg	3,6 kg
Watt	50 W	50 W	26 W
Operation	Touch control	Touch control	Touch control
Brush/polish function	yes	yes	no
Max. speed	1000 rpm	1000 rpm	⊗
Adjustable speed	20 positions	20 positions	⊗
Speed control	Touch control	Touch control	⊗
Brush 20 mm soft	● 51403	● 51403	⊗
Brush 45 mm soft	● 51409	● 51409	⊗
Brush 45 mm hard	● 51408	● 51408	⊗
Pumice stone 45 mm	● 51405	● 51405	⊗
Vacuum function	no	yes	Nee
Adjustable vacuum	no	yes	Nee
Glass capillary for comedones	⊗	● 51421	⊗
Glass capillary for wrinkle correction	⊗	● 51422	⊗
Glass capillary for small beaker 1 cm	⊗	● 51423	⊗
Glass venole large beaker 3 cm	⊗	○ 51424	⊗
Glass venole large beaker 4 cm	⊗	○ 51428	⊗
Spray function	no	yes	no
Iontophoresis function	no	no	yes
Adjustable voltage	⊗	⊗	20 positions
Ionisation	⊗	⊗	yes
Desincrustation	⊗	⊗	yes
Spherical/roll applicator	⊗	⊗	● 61461



High-frequency machine

The Bentlon high-frequency machine has an activating effect on the skin. The skin is disinfected and circulation is stimulated. This makes it particularly suitable for use on greasy skin, blocked skin and acne. The various glass electrodes allow treatment of face and body, whereby the desired strength can be set and read on the digital display.

Cooling/heating unit


This compact high-tech cooling/heating unit can be used for practically every treatment. The 14-carat gold-plated side means that it can be heated from 32 to 42 degrees. The heat can be used to stimulate, regenerate and vitalise the skin. The cell metabolism is also stimulated, more oxygen is transported to the cells and the fibroblasts are stimulated to produce collagen. Active ingredients in creams, masks, ampoules or concentrates are also absorbed better by the skin. The silver-coloured inox steel side can cool to 5 till 10 degrees, so that the skin is sealed after the treatment. It provides a cooling and fresh feeling, calms the skin after an intensive treatment and reduces any swelling or pain.




Ref. 51485




Ref. 51480

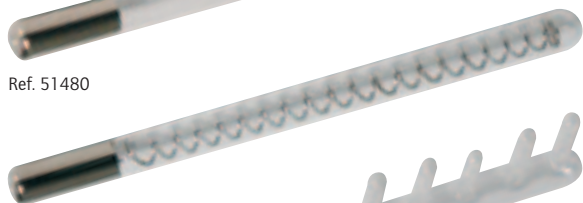

Ref. 51486

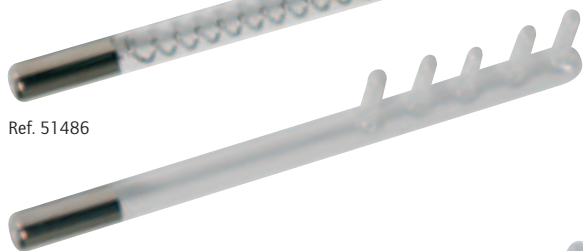

Ref. 51484

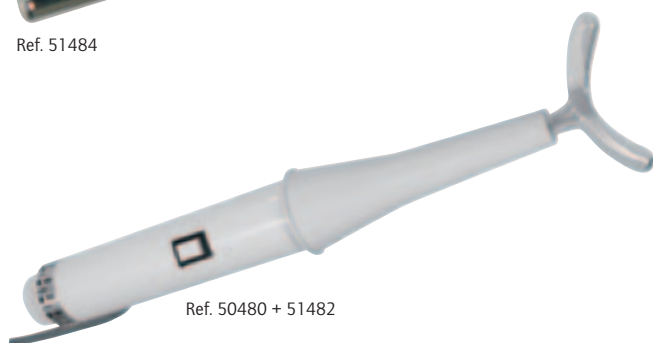

Ref. 50480 + 51482



Ref. 60554

Device	High-frequency handpiece
Ref. 230 V	50480
Ref. 115 V	⊗
Width x depth x height	320 x 45 mm
Weight	1,5 kg
Watt	30 W
Operation	Turning control on holder
Adjustable voltage	10 positions
Face electrode 20 mm	○
Face/body electrode 45 mm	○
Neck electrode (curved)	○
Comb electrode	○
Rod/ spiral electrode	○

Device	Thermostick
On basic version	○
Ref.no.	60554
Voltage	230 V
Consumption	40 W
Cooling	5 - 10 degrees
Cooling adjustable	5 stages
Heat	32 - 42 degrees
Head adjustable	5 trappen
Display	LED

Accessories	Ref.
Kam-elektrode	51484
Rod / spiral electrode	51486
Rod / spiral electrode 20 mm	51480
Rod / spiral electrode neck	51482
Rod / spiral electrode 45 mm	51485

BENTLON



Ref. 50560



Ref. 50580



Ref. 50570

Ultrasonic cleaner Silver

The HAM Code lays down guidelines for optimum cleaning and disinfection. The Bentlon Ultrasonic cleaner Silver complies with these guidelines. The cleaning time can be set to 90, 180, 280, 380 or 480 seconds and the digital display indicates the remaining cleaning time.

Ultrasonic cleaner Gold

The Bentlon Ultrasonic cleaner Gold has a bigger capacity than the Silver and this means that it is possible to clean larger instruments and more instruments at the same time. In addition to the same time settings as the Silver, the Gold also offers the option of heating the liquid to 60 degrees so that stubborn or dried on dirt is removed.

Ultrasonic cleaner Platinum

The Bentlon Ultrasonic cleaner Platinum is the perfect product to professionally clean and disinfect your instruments. The Platinum version has a 3 color LED display which enables the user to check the temperature and remaining cleaning time. The cleaning time can be set up to 30 minutes at 5 different temperature levels. Moreover, the device is protected against overheating and switches off automatically after 90 minutes. The Bentlon Ultrasonic cleaner Platinum has an extended tank capacity of 3 ltrs and allows to place products of up to 26,6 cm. long. The device is practical and safe in use.

Hot Cabinet Bronze

The Bentlon Hot Cabinet Bronze steams and heats towels to the correct temperature. The towels stay hot for up to 60 minutes with the timer. The Hot Cabinet Bronze also has a useful read-out function for the water reservoir. This means that you can see exactly how much water to add to heat the towels.



Ref. 50620

Hot Cabinet Silver

The Hot Cabinet Silver holds more than the Bronze so that you can heat more towels at the same time. The Silver offers the facility to heat damp towels directly. This means there is no water reservoir.



Ref. 50625

Towel heater	Hot Cabinet Bronze	Hot Cabinet Silver
Ref.no.	50620	50625
Function	Towel heater	Towel heater
Voltage	230 V	230 V
Watt	850 W	200 W
Width x depth x height	21 x 23,3 x 31,5 cm	45 x 31 x 36 cm
Capacity	± 8 pcs.	± 50 pcs.
Heating time	± 10 min	± 15 min
Timer	●	⊗
Automatic stop time	●	●

Desinfection/sterillisation equipment	Ultrason cleaner Silver	Ultrason cleaner Gold	Ultrason cleaner Platinum
Ref.nr. 230 V	50560	50570	50580
Function	Desinfection	Desinfection	Desinfection
System	Ultrason	Ultrason	Ultrason
Width x depth x height	23 cm /18cm /16 cm	26,4 x 16,4 x 8 cm	25 x 17,4 x 7,8 cm
Watt	70 W	170 W	170 W
Frequency	●	●	●
Desinfecting bowl	⊗	⊗	⊗
Timer	●	●	●
Temperature indicator	●	●	●
Pressure	⊗	⊗	⊗
Automatic stop time	●	●	●
Automatic fault alert	⊗	⊗	⊗
Water drainage	⊗	⊗	●

BENTLON

Wellness



Touch of Excellence



Wellness couches	48
Wellness equipment	52



BENTLON

gold

Wellness couch



Ref. 60040

Comfort, ease of use and countless options

- Luxurious design
- Optimal comfort
- Available in various designs
- With option for heating and massage

It is with pride that Bentlon presents to you a series of Wellness Couches. The Wellness Couch features an exceptionally high comfort level and luxurious upholstery and frame finish. The Wellness Couch is not only ideal for giving body treatments but can be just as easily used for facial treatments. More comfort, user-friendliness, style and countless possibilities - get this and more with the Wellness Couch!



The Wellness Couch Gold can be adjusted to a very low position (60 cm) to allow easy access for clients. The couch can then be raised electronically to the ideal working height (up to 90 cm). This allows both short and tall specialists to adopt a good work posture when providing treatments. The back and leg rests can also be adjusted electronically. This makes it possible for both you and your client to maintain the most comfortable position during the entire treatment. An electric foot control ensures smooth and easy adjustment. The Wellness Couch is also equipped with four small wheels for easy handling.



Special 10 cm-thick high-quality elastic foam filling makes the couch extremely comfortable. The upholstery consists of seven sections, including adjustable armrests and a patented Bentlon head support that can be used when the client is lying flat.



In addition to the standard options offered by the Wellness Couch, there are a number of additional options that can be selected. One such option is seat heating. This is not only especially pleasant for clients, but also energy saving because the temperature in the salon can be set a few degrees lower, creating a pleasant working climate for the specialist. Another option is the built-in massage system. This system features three different programmes that vary in intensity, and with the additional option of partial massages. A must-have for every couch!

Arm support

As the rest of the couch, the arm support has a 10 cm-thick filling and can be adjusted continuously up to 14 centimetres.











The Bentlon design and style of this new Wellness Couch is truly amazingly beautiful. The casing is constructed of solid wood and finished with a high-quality HPL top layer. This makes it not only beautiful to look at, but hygienic, easy to clean and both scratch and colour fade-resistant. The casing is also designed to prevent any 'views' of the underlying frame.

Choose your favourite style

The wooden casing is available in six different colours, varying from white and light wood, to medium-dark, dark and black wood. The couch upholstery is available in white, black, Dakar beige and grey and is made of high-quality polyurethane imitation leather. The surface is both easy to clean and waterrepellent. It retains its beauty, suppleness and shape, even when adjusted back to the starting position. You can put together your own unique Wellness Couch by combining your favourite colours for the wooden casing and upholstery. All combinations are possible!

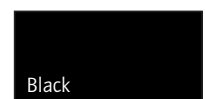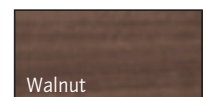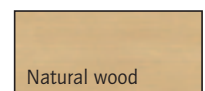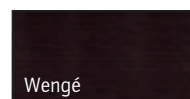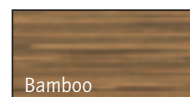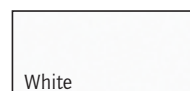

Ergonomically designed for specialists

The new Bentlon Wellness Couch is specially designed to offer exceptional comfort to both you and your clients. Because the couch is height adjustable, you have plenty of leg space when doing facial treatments. You can also comfortably rest your arms on the head support when giving treatments. This means you no longer have to bend over and can work standing upright, promoting a more comfortable and relaxed work posture. When giving body treatments, the specially designed wooden casing gives you extra foot room, allowing you to sit closer to your client and, once again, reducing pressure on your back.

Portable Carrying table

The massage table for ambulant work. The new Bentlon portable carrying table only weighs 16 kg and is easy to fold down and out from bag to table. The table is not only light, but also looks good. The table has a head rest that is easy to attach.

Wellness Couch Bronze

When the massage table is used intensively, electrical adjustment of the back support is desirable, particularly for older or heavier clients. The Wellness Couch Bronze massage table has an electrical adjustment of the back support that is continuously adjustable by means of a hygienic foot operation.



Ref. 50106.00



Ref. 60115.100



Ref. 60040.100.02

Wellness couches	Portable carrying table	Wellness Couch Bronze	Wellness Couch Gold
Ref.no.	50106	60115 (230 V) 60116 with heating (230 V) 60117 (115 V)	60040 60041 (with heating) 60042 (with heating and vibra)
Length	184 cm	200 cm	220 cm
Width	70 cm	67 cm	84 cm
height	60 - 81 cm	70 - 130 cm	60 - 94 cm
Weight	16 kg	80 kg	150 kg
Max. permitted weight	150 kg	150 kg	180 kg
Upholstery	Imitation	Imitation	Imitation
Coloured upholstery	White	37 colours	4 colours
Back adjustment	⊗	●	●
Height adjustment	Manual	Electrical	Electrical
Height control	Manual	Foot switch	Hand control
Leg support adjustment	⊗	⊗	Electrical
Breathing gap	●	●	●
Vibra	⊗	⊗	○
Chair heating	⊗	○	○
Case material	⊗	⊗	Wood
Wood models	⊗	⊗	4 wood/2 colour
Foldable arm support	⊗	⊗	●
Memory reclining and rising pos.	⊗	⊗	●

For all the colors of the upholstery go to page 76

BENTLON

Wellness equipment



To make your salon complete!

- Offers customers ultimate relaxation
- Professional hotstone equipment

Hotstone heater Silver

The Bentlon Hotstone Heater Silver is a smaller version of the Gold and has a capacity of 5.5 litres. Just like the Gold, the Silver also has automatic temperature control and the heat can be read off the separate thermometer.

Hotstone heater Gold

The hotstone massage relaxes the mind and body. The Bentlon Hotstone Heater Gold has a large capacity for heating the stones. The heater has programmed temperature settings and automatic temperature control. The temperature can be read off the separate thermometer.

Basalt stones

The basalt stones have been polished by hand so that all the stones are the same shape and there are no pressure point differences. There are two kinds of basalt stones. There is a set with 36 stones and a set with 64 stones. Of course, the stones are different in size so that you can massage any body part required. The stones are supplied in a luxury wooden storage box.



Ref. 60625



Ref. 60630



Ref. 61025

Hotstone heater	Silver	Gold
Ref.	60625	60630
Length	32 cm	50 cm
Depth	22 cm	35,5 cm
height	10 cm	32 cm
Capacity	5,5 liter	20,5 liter
Heating-up time	15 Minutes	15 - 20 minutes
Voltage	230 V	230 V
Power consumption	200 W	550 W
Thermostat	Yes	Yes

Basalt stones	36 pcs	64 pcs
Ref.	61025	61026
Extra large stones	⊗	2
Large stones	8	14
Middle stones	12	14
Small stones	16	7
Trigger stones	⊗	4
Toes	⊗	8
Chakra stones	⊗	7
Oval stones (eyes)	⊗	2
Small oval stones white marble	⊗	6



Ref. 60500



Ref. 60505



Ref. 60510

Completely in line with the developments of the Wellness branch, Bentlon is launching heat blankets for body packs, the Bentlon Heat blanket. In a body treatment the pack, depending on the kind used, should be left to work in for about 20 minutes. During this time the client should enjoy a luxurious sensation from the pleasant warmth. Just covering with a blanket or duvet is not sufficient. With the Bentlon heat blanket you create this luxurious heat. Moreover, the substances of the pack are absorbed better by the skin due to the heat. The blankets are made of waterproof plastic. They are easy to clean and are water tight which means that they are very hygienic to use.

Heat Underblanket

The Bentlon Heat Underblanket is an underblanket that the client lies on. It can be attached to the treatment chair or table and has a comfortable head cushion. The underblanket can be set to 3 temperature levels and has a digital timer that gives a signal at the end of the treatment time. The Bentlon Heat Underblanket is also highly suitable for use as comfortable chair heating during normal treatments.

Heat Overblanket

The Bentlon Heat Overblanket is a heat blanket that you place over the client. It has the same temperature control and timer as the underblanket described above. The overblanket is ideal to replace a duvet. Your client is nice and warm and you can set the temperature in the salon lower which will give you a more comfortable working climate.

Heat Coverblanket

The Bentlon Heat Coverblanket is an extra large blanket the client lies on and is wrapped in. The cover blanket has three heat zones that can be set separately, i.e. upper body, lower body/upper legs and lower legs. Completely according to the type of treatment or pack you can set the temperature of each zone individually.

(All Bentlon Heat blankets have a patented Electro Block safety system and operate on 24V low current.)

Heatwave

The infrared beams of the Heatwave lamp deeply penetrate the skin. They stimulate blood circulation and the absorption of skin care products is intensified. The Heatwave also has a very relaxing effect. Infrared light is also used on a large scale to relieve pain and to treat infections, skin damage, muscle aches and cellulite.

The Heatwave support is constructed in such a way that it can be placed under a massage table, but it also fits around the stand of a treatment chair. It is adjustable in height and the lamp can also be tilted.

There are two versions of the Heatwave, Heatwave Silver and Heatwave Gold. In the Heatwave Silver version the infrared fitting can be attached to the ceiling and the height can be adjusted easily. The Heatwave Gold is supplied on a stand.

Heatwave	Gold
Ref.no.	60495
Light unit dimensions	100 x 32 cm
IR lamp reflector dimensions	72,5 x 12 cm
Stand Dimensions	102 x 92 cm
Height	120 – 160 cm
Weight	28 kg
Watt	1300 W
Type of light	IR
Frequency	600 – 1200 nm
Treatment surface	80 x 100 cm
Mobility	4 castors
Timer	30 minutes
Automatic stop time	●



Ref. 60495

Heat blankets	Underblanket	Overblanket	Coverblanket
Ref.no.	60500	60505	60510
Dimensions	160 - 170 cm	160 - 170 cm	160 - 170 cm
Weight	6,5 kg	6,8 kg	13,5 kg
Watt	253 W	253 W	253 W
Type	Under blanket	Over blanket	Packable over blanket
Material	Plastic	Plastic	Plastic
Power supply	24 V	24 V	24 V
Safety	Electro Block	Electro Block	Electro Block
Watertight	●	●	●
Timer	●	●	●
Automatic stop time	●	●	●
"Time up" alert	●	●	●
Heat adjustment	Whole blanket	Whole blanket	3 sections

BENTLON

Podo



Touch of Excellence



Podiatry chairs	58
Podiatry devices	66
Podiatry S Cabinet	70
Cryopen by Bently	72
Podiatry accessories	74



BENTLON



Podo



Ref. 65814.120

...without compromising Bentlon quality!

- Functional design at entry-level pricing
- Optimal ease of use because of 2 actuators
- High working position
- Lifting capacity of 200 kg
- Innovative one touch treatment positions
- Ergonomic high quality upholstery





The Bentlon Bronze treatment chairs are equipped with 2 actuators which make it possible to electrically adjust the height and tilting function of the chair to set it to the desired position. It is also possible to put the chair in massage position electrically. The back rest adjustment is easy to operate by means of a gas compression spring. This combination of adjustment possibilities offers the specialist exceptional ease of use.



The treatment position of the Bentlon Bronze ensures a comfortable seated position for the client at all times. Thanks to the narrow back rest, the chair has a slender appearance, which makes the salon or consultation space appear more spacious. In addition, the arm supports can fold up completely, which makes it easy for the client to step into and out of the chair from the side.

To give the Bentlon Bronze chairs a luxurious appearance, Bentlon has selected high-quality imitation leather with a calfskin grain and a subtle matte-lustre finish. The upholstery is very easy to clean and is resistant to disinfectants. The abrasion- and tear-resistant upholstery ensures a longer service life.



Clients can step into and out of the chair easily by using the unique, patented electric step-in and out seat assist. This step-in and out seat assist lowers the chair, tilts the seat forward and rotates the leg and foot support inward. The chair follows and supports the client's movement as they step in or out. The chair also features robust construction which significantly increases the lifting and tilting power.

Extra-low step-in seat assist

The extra-low step-in seat assist makes it even easier to step into and out of the chair. The seat of the chair can be lowered to a height of 51 cm. For elderly or people who have trouble with mobility, this provides good support—only a little effort is required to step into and out of the chair.

High working position

The Bentlon Beauty Bronze is characterized by an exceptionally high working position. The seat of the chair can in fact be raised to a height of 85 cm, which also allows the specialist to give a massage at this working height. If the specialist wants to perform a foot treatment, then that is also possible. In that case, the working height can be as high as a maximum of 93 cm. This allows the specialist to comfortably work on the client's feet while standing.

Foot switch

It is also possible to operate the Bentlon Beauty Bronze using a foot switch. This makes it easy to change position during the treatment, while maintaining hygienic working conditions.

BENTLON



Podo Rotate



Ref. 65880.123.00

180° rotatable podiatry chair

- **User-friendly because of its ability to rotate**
 - Ideal when getting in or out of the chair
 - For maximal use of practice space
- **Electrical height, tilting and back adjustment using a foot control**
- **Working height of max 130 cm**
- **High quality upholstery for excellent seating comfort**
- **Ergonomically shaped foot rests**
 - Foldable foot rests
 - Independently adjustable

The Bentlon Podo Gold Rotate is the ideal solution for smaller practices. Being able to rotate this professional podiatry chair makes it the perfect match for podiatrists looking to maximally use practice space whilst customers can still easily get in and out of the chair. The combination of foldable arm rests and the option to rotate allow patients to get in the chair without too much effort.



Additionally, much attention is paid to the ease of use for the podiatrist. By using a foot control height, back and tilting position can be adjusted, reaching a maximum working height of 140 cm. The chair can easily be rotated to the left or right using the lever below the seat.



In terms of ergonomics the Bentlon Podo Gold Rotate meets the highest quality standards. The chair is covered with luxurious leather upholstery available in 37 colors and contains a special foam for optimal seating comfort, up until the foot rests. The footrests moreover, are foldable which is very practical when using a foot bath. The footrests can also be adjusted independently from each other.



Ref. 65880.123.00

BENTLON



Podo



Ref. 65864.120.02

Reaching revolutionary heights!

- Revolutionary working height of 153 cm
- Optimal ease of use because of 3 or 4 actuators
- Patented Step Out - Assist
- Home position and memory position
- Lifting capacity of 200-225 kg
- Ergonomic upholstery available in 37 colors!



The Bentlon Podo Platinum is the most advanced treatment chair in its class! This pedicure chair combines all of the aspects needed by a professional pedicurist to be able to work efficiently and comfortably. Very high working height, stability, ease of use and extremely high quality!

Working height

Probably the most important item in a professional pedicure practice is the treatment chair. Since everything revolves around the treatment of the feet, it is an absolute must to have the client's feet at the ideal height, for the client to sit comfortably while the pedicurist can assume an ergonomic working posture. The Bentlon Podo Platinum treatment chair fulfills this requirement and offers a working height of up to 144 cm.

Additional working height can be reached with the models equipped with an extendable frame. With this innovative solution a working height of up to even 153 cm is possible! That makes this treatment chair 25% higher than its predecessor!



Optimal ease of use for the pedicurist has been the most important starting point in the development process. Therefore this treatment chair is equipped with a minimum of 3 motors providing electric adjustments for the height, back and seating angle of the chair. Using a foot control, the treatment chair can be adjusted from the home position to the desired treatment position and then, with just one press of the pedal, back again to the home position.

Furthermore, the foot control has an extra function for setting a memorized, preferred treatment position, which then can be initiated with just one press of the pedal. This way it is for example possible to automatically set the chair to the Trendelenburg position. Software settings ensure that the chair can never be placed into unwanted positions.



Ref. 65860.124.02



The Bentlon Podo Platinum is equipped with a patented step-in and step-out seat assist that supports the client while stepping in or out of the chair. This seat assist raises the chair and tilts the seat cushion forward. The chair follows and supports the client's movement as they step in or out. This step-in and out seat assist can also be activated with the foot control.



The Bentlon Podo Platinum offers optimal comfort because its upholstery is ergonomically contoured according to the body shape which provides the client with a very comfortable seated position throughout the entire treatment. The upholstery is made of high-quality imitation leather with a calfskin grain resulting in a very classy, luxurious look. Moreover, the Bentlon Podo Platinum can be customized to everyone's desire as it is available in 37 different colours upholstery and 2 colour choices for the foot frame.

Stability

The Podo Platinum is the most stable treatment chair ever made and can lift weights of up to 210 - 225 kg. This enormous load bearing capacity means that the chair is suitable for a wide range of clients. The frame can also be anchored to the floor or can be fitted with wheels making the chair portable.

Lower Arm Support



Together with the Delft University of Technology, Bentlon has made an important invention with a worldwide patent: the Lower Arm Support. This is an ergonomic support with bin for podiatrists and represents a breakthrough in the podiatry industry. The work of the podiatrist requires a high level of precision, concentration and long periods of muscle contraction. During the treatment a certain posture is adopted or the load on the body is on one side. Often podiatrists sit with a bent posture or in a slightly twisted position to obtain a good view when working and their arms are continuously raised while working. This is very tiring and wearing for the body.

Support

The Lower Arm Support supports the lower arm so that the weight of the upper arm, lower arm and hand does not have to be supported wholly by the muscles in the back and shoulders. Because of this support the treatment chair that the patient sits on can be set higher; this means that the podiatrist can see the foot better and is able to adopt a more upright position of the lower back and the neck. This means a considerable improvement in the work posture and is much less tiring for the shoulders and arms.

Hygienic

The Lower Arm Support also has a removable bin that is easy to remove and put back and is easy to clean due to its funnel shape. The support bracket and the bin have been developed so that there is enough space around the foot to work on it. The Lower Arm Support is suitable for left-handed and right-handed podiatrists.



Ref. 65814



Ref. 65880




Ref. 65864

Podo chairs	Bronze	Gold Rotate	Platinum	Platinum extendable bridge
Ref.no. with hand control	65814	⊗	⊗	⊗
Ref.no. with foot control	65817	65880	65860	65864
Ref.no. 1 leg support	⊗	⊗	65862	65866
Number of actuators	2	3	3	4
Length	170 cm	178-200	197-199 cm	197-199 cm
Chair width	75 cm	80 cm	78 cm	78 cm
Sitting width	55 cm	61 cm	53 cm	53 cm
Height	51-83 cm	52-107	53-106 cm	53-115 cm
Chair weight	53 kg	74 kg	65 kg	78kg
Max. permitted weight	200 kg	200 kg	225 kg	210 kg
Max. working height	93 cm	130 cm	136 cm	145 cm
Watt	96 W	96 W	85 W	85 W
Upholstery	Imitation leather	Imitation leather	Imitation leather	Imitation leather
Coloured upholstery	3 colours	37 colours	37 colours	37 colours
Case colour white/grey	●	○	○	○
Extendable head support	●	⊗	⊗	⊗
Tilting head support	⊗	⊗	⊗	⊗
Head support depth adjustment	⊗	⊗	⊗	⊗
Extendable leg support	●	●	●	●
Fold-away arm supports	95	185	185	185
Fold-away leg supports	⊗	●	●	●
Lying and exit position	●	⊗	●	●
Memory reclining and rising pos.	⊗	⊗	●	●
Rotatable	⊗	●	⊗	⊗
Height adjustment	Electrical	Electrical	Electrical	Electrical
Tilt adjustment	Electrical	Electrical	Electrical	Electrical
Back adjustment	Gasdrukveer	Electrical	Electrical	Electrical
Leg support adjustment	Pull-override	Pull-override	Pull-override	Pull-override
	Gas compression spring	Gas compression spring	Gas compression spring	Gas compression spring

For all the colors of the upholstery go to page 76


Accessories	Ref.
Chair covering towelling	61849
Protective set head/arm/foot	61859
Lower Arm Support narrow	66068/61854
Lower Arm Support wide	66069/61088/61853
Support roll round	61090
Support roll semi-circular	61091
Leg support cushion	66067
Castors	61099
Upholstered cushions	66066
Massage head support	61083
Cover massage head support	61022
Chair cover	65860/65864
Chair cover	65862/65866/61860/61862

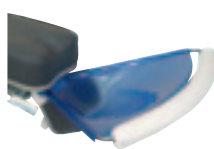

Ref. 61088



Ref. 61853

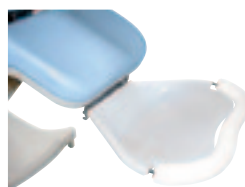

Ref. 66068



Ref. 61854

BENTLON

Podiatry devices



Experience a powerful motor and a lightweight hand piece!

- Innovative models for suction or spray technique
- Powerful motors of up to 40.000 RPM
- Light weight hand pieces that prevent fatigue
- Awarded the Dutch award for Good Industrial Design



Bentlon launched a new line of podiatry devices, designed specifically to meet the needs and requirements of professional pedicurists and podiatrists. The line comprises 7 spray and 2 suction models, each of which is fitted with an ergonomic, lightweight hand piece and designed with hygiene and quality as the starting point. The new design of the podiatry devices will be a real eye-catcher in your practice!

Innovative spray technology

The podiatry devices with spray technology are silent motors with a powerful range of between 5,000 and 30,000-40,000 RPM. Up to 14 or 19 speed settings can be chosen within this range, depending on the type of machine, and the unique spray technology ensures that both the skin surface and the bit stay cool and that a fine mist is sprayed. All of the models have a large, transparent 245ml reservoir which makes it very easy to keep an eye on the level of disinfectant spray.

Proven suction technology

The suction-powered podiatry devices are fitted with very powerful motors that have a maximum suction power of 190 litres per minute and have the same range as the spray models. The speed can be varied with 14 or 19 settings and the dust bag doesn't need to be replaced after every use. When required, an indicator is shown via an automatic timer on the LED display.

Lightweight hand pieces

Every podiatry device is fitted with a professional hand piece from the Smart or Perfect series, both of which have unique properties. All of the hand pieces are remarkably light and therefore prevent fatigue and strain.

These ultra-light hand pieces have been designed with hygiene, sterilisation options, ergonomics and user-friendliness in mind and the result are hand pieces that maximise comfort and are easy to use.

Spray Smart B
Ref. 65150



Spray Smart
Ref. 65155



Spray Smart BL
Ref. 65160



Spray Perfect /L
Ref. 65170



Spray Perfect BL
Ref. 65195



Vacuum Smart B
Ref. 65100



Vacuum Smart B
Ref. 65105



The hand pieces from the Smart series are a single unit with an integrated spray or vacuum bit head with micro-motor. The Smart hand pieces for the vacuum models are fitted with 'soft touch pads' for an optimum and comfortable ergonomics. Vibrations are absorbed by the innovative way in which the micro-motor has been integrated into the hand piece. A knob on the models with a Smart series hand piece can be used to adjust the quantity of spray or suction power and touch keys can be used to adjust the RPM.

The hand pieces from the Perfect series comprise a separable spray bit head and micro-motor, making it possible to clean the bit head separately and sterilise it in an autoclave. It is also enables to fit a Bentlon Contra Angle or Bentlon Omnicut onto the hand piece, so that the hand piece can be used for multiple purposes. These hand pieces are made of titanium, making the surface very scratch-resistant and comfortable to hold. Both the RPM and the quantity of spray can be set using the touch key panel.

Both the Smart and Perfect series include hand pieces with LEDs, which provide extra lighting at the treatment site and give the pedicurist or podiatrist a better view of the cause being treated.



Ref. 65195

Ref. 65105

Ease of use

In addition to the compact design, ultra-light weight of the hand pieces, ease of cleaning and different options for every model, the podiatry devices with a touch key panel offer the option to configure three preferred spray/vacuum settings and three preferred speeds. This enables the pedicurist or podiatrist to configure the programmes to his/her personal method of working and is a very practical functionality for starting the treatment quickly. All of the models also have a handle that practically retracts into the machine when not folded out and is an ideal extra for ambulant work. The models that are fitted with a brushless motor also require minimum maintenance.

Exclusive appearance

In addition to the basic functionality, the sophisticated models have extra features such as LEDs in the machine, which give it a very exclusive look and contribute to the professional appearance of the podiatry device.

Accessories

Extra accessories are available for every model. All of the models can be supplied with a foot pedal, bits and panaspray.

Dutch design

Bentlon designs all of its products according to the latest technologies and combines modern design with high quality and reliability. The new podiatry devices are fitted with high-tech motors and are wonderful trendy eye-catchers, thanks to their clean, modern design and high-gloss finish.

And the digital display, touch key panel and LEDs give the machine an extra touch of prestige.

Quality and service

The new podiatry devices meet all of Bentlon's high quality requirements. They are fire-proof and fitted with high-quality motors that have a very long service life. Bentlon has been a professional manufacturer of pedicure and beauty care equipment for almost 40 years and as such is market leader within its sector. Bentlon stands for high quality and provides a full two years' warranty on its products, and even five years on manufacturing defects. Moreover, it is possible to take out a service contract on the pedicure machines, in order to guarantee the life service of your product for even longer.

equipment	Smart B	Smart	Smart BL	Perfect	Perfect	Perfect L	Perfect BL
Ref.no.	65150	65155	65160	65170	65180	65190	65195
Voltage	autorange	autorange	autorange	autorange	autorange	autorange	autorange
							
Hand piece	Smart B	Smart	Smart BL	Perfect	Perfect	Perfect L	Perfect BL
Weight hand piece	112 gr	105 gr	98 gr	148 gr	148 gr	151 gr	129 gr
Brushless motor	Yes	No	Yes	No	No	No	Yes
Power	120W	120w	120W	120W	120W	120W	120W
LED	No	No	Yes	No	No	Yes	Yes
Speed	30.000 rpm	40.000 rpm	40.000 rpm	40.000 rpm	40.000 rpm	40.000 rpm	40.000 rpm
Speed settings	14	19	19	19	19	19	19
Dimensions (WxDxH) cm	30,2x22,5x17	30,2x22,5x17	30,2x22,5x17	30,2x22,5x17	30,2x22,5x17	30,2x22,5x17	30,2x22,5x17
Weight	3,4kg	3,4kg	3,4kg	3,4kg	3,4kg	3,4kg	3,4kg
Display	LCD Blue	LCD Blue	LCD Blue	LCD Blue	LCD Blue	LCD Blue	LCD Blue
Spray reservoir (ml)	245 ml	245 ml	245 ml	245 ml	245 ml	245 ml	245 ml
Reservoir	Removable	Removable	Removable	Removable	Removable	Removable	Removable
Operated by	Turnable knob	Turnable knob	Turnable knob	Turnable knob	Touch control	Touch control	Touch control
Spray settings	Variable	Variable	Variable	Variable	10	10	10
Memory spray	n/a	n/a	n/a	n/a	3	3	3
Memory speed	3	3	3	3	3	3	3
Noise (db)	50dB	50dB	50dB	58dB	58dB	58dB	58dB
Carry handle	Yes	Yes	Yes	Yes	Yes	Yes	Yes
Footcontrol option	Yes	Yes	Yes	Yes	Yes	Yes	Yes
LED Bentlon logo	No	No	No	No	Yes	Yes	Yes

Vacuüm equipment	Bronze Smart B	Bronze Silent Smart B	Silver Smart B	Silver Silent Smart B	Accessories	Ref.
Ref.no.	65100	65110	65105	65112	Foot switch with cable	60856
Voltage	230V *	115-230V	230V	115-230V	PanaSpray (oil for hand piece)	66070
					Dustbag for vacuum models	66078
Hand piece	Smart B	Smart B	Smart B	Smart B	Wireless foot switch	66065
Weight hand piece	112 gr	112 gr	110 gr	110 gr	Contra Angle	66075
Brushless motor	Yes	Yes	Yes	Yes	Contra Angle with light	66076
Power	230W	90W	230W	90W	Omnicut (to remove chapped skin)	66060
LED	No	No	No	No	Blades for Omnicut	66061
Speed	30.000rpm	30.000rpm	40.000rpm	40.000rpm		
Speed settings	14	14	19	19		
Dimensions (WxDxH) cm	30,2x22,5x17	30,2x22,5x17	302x225x170	302x225x170		
Weight	3,4kg	3,2kg	3,4kg	3,2kg		
Display	LCD Blue	LCD Blue	LCD Blue	LCD Blue		
Reservoir	Dust bag	Dust bag	Dust bag	Dust bag		
Reservoir	Removable	Removable	Removable	Removable		
Operated by	Touch control	Touch control	Touchcontrol	Touchcontrol		
Max. suction power	190 ltr/min	200 ltr/min	190 ltr/min	200 ltr/min		
Memory vacuum	3	3	3	3		
Memory motor	3	3	3	3		
Noise (db)	65dB	48dB	65dB	48dB		
Carry handle	Yes	Yes	Yes	Yes		
Footcontrol option	Yes	Yes	Yes	Yes		
LED Bentlon logo	No	No	No	Yes		
* Ref.no. 115V	65102		65107			



Ref. 66060

Podo S Cabinet



The ultimate working station for larger practices

- Completely adjustable to your wishes
- Various designs
- Choice between different hand pieces

Podo S Gold Cabinet

The Podo S Gold Cabinet has two spacious storage drawers and a drawer with UVC lighting to keep instruments sterile. A magnifying lamp holder can be positioned on the worktop in four different places. The motorised arm and motor unit itself can be rotated so that they can be moved to the most convenient working position. After use, the arm can be folded back under the worktop. The Podo S Gold Cabinet Base ranks among the best equipment available in the professional pedicure industry.

Bentlon Podo S Platinum Cabinet

The Bentlon Podo S Platinum is the ultimate work station for larger pedicure practices, offering numerous possibilities. It is not a standard pedicure unit, but is built according to your specific requirements. If desired, the Podo S Platinum can be fitted with the Pododent 1 or Pododent 2 hand attachment with or without light. The filing head of the Pododent 2 hand attachment is removable. The hand attachment can then be converted into a contra-angle with a speed of 200,000 revs using a special accessory, or into an electrical hard skin file using the Omnicut attachment. The Podo S Platinum can also be fitted with an air-driven turbine hand attachment with a speed of 430,000 rpm,

with or without light as required. The unit comes with a separate spray/compressed air hand attachment as a standard accessory. Extraction is also available as an option. The Podo S Platinum can be integrated into a luxurious mobile cabinet.

Compressors

The Podo S Platinum must be equipped with a compressor in order to use the spray function the hand attachments and the spray and compressed air functions on the syringe. If you choose to have the Podo S Platinum integrated into a cabinet, it will be fitted with a built-in silent compressor.

Podo equipment Cabinet

Options (available in 7 colours)

Ref.

Bentlon Gold Cabinet	65625
Bentlon Gold Cabinet + UV drawer	65626
Bentlon Gold Cabinet with arm + tableau	65627
Bentlon Gold Cabinet met arm + tableau + UV lade	65628
Bentlon Platinum Cabinet M + S	65631
Bentlon Platinum Cabinet M + S + UV drawer	65636
Bentlon Platinum Cabinet M + S + Turbine	65632
Bentlon Platinum Cabinet M + S + Turbine + UV drawer	65637
Bentlon Platinum Cabinet M + S + vacuum	65633
Bentlon Platinum Cabinet M + S + vacuum + UV drawer	65638
Bentlon Platinum Cabinet M + S + vacuum + Turbine	65634
Bentlon Platinum Cabinet M + S + vacuum + Turbine + UV drawer	65639

Choose your handpiece/motor (only on Platinum Cabinet)

Pododent 1 motor	65710
Pododent 2 motor without light	65711

Choice of turbine (only on Platinum Cabinet)

Turbine motor Light (400.000 rpm) Bentlon podo Platinum	65722
---	-------

Compressors (included in Platinum Cabinet as standard)

Air	32 l/min
Pressure	8 bar
Noise level	69 dB


Accessoires

Ref.

Bentlon Podo Omnicut (only suitable for Podent 2 hand attachments)	66060
Bentlon Podo Omnicut blades	66061
Draadloze foot control for Podo Gold + Platinum	66065
Contra-angle (only suitable for Podent 2 hand attachments)	66075
Contra-angle met licht (only suitable for Podent 2L hand attachments)	66076
Aluminium tray anodised	66071
Aluminium instrument rack	66072
Plastic disposable tray	66073
Panaspray	66070



Ref. 65625.06

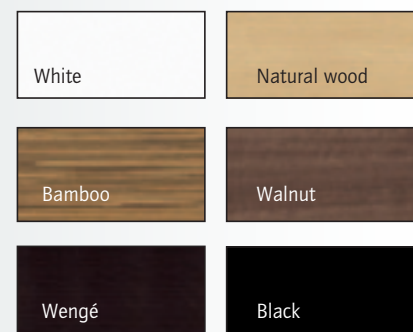

BENTLON

Cryopen by Bentlyon



Ref. 65300

Removes warts quickly,
effective and virtually painless!

- Quick procedure: Only 1 – 2 minutes!
- Sterile procedure: No cutting, no bleeding, no infection risk
- Cold temperature reduces skin sensibility almost completely
- Painless treatment
- No unwanted side effects for the patient
- Safe in use, no scars
- Excellent results

Removing foot warts is very precise, effective and almost painless with the CryoPen by Bentlon. The very fine beam of liquid gas of the CryoPen by Bentlon ensures that the treatment can be carried out precisely by millimeter leaving the healthy skin around the wart untouched. Cryotherapy is already a popular technology used by many doctors, dermatologists and pedicurist.

Cryotherapy

Cryotherapy is the latest technology in the treatment of warts. Because of controlled freezing warts can be removed in a very quick procedure of only 1-2 minutes.

Cryotherapy has many benefits compared to traditional treatments. First of all, a treatment with the Cryopen is so precise that healthy skin tissue around the wart remains intact. It also rules out blistering. Moreover, skin sensitivity is minimized because of the cold temperature which makes the treatment almost completely painless. Because of the constant temperature of -89°C and the pressure of the CryoPen treatments can be executed very precisely and effectively.

The procedure using the CryoPen by Bentlon is also very sterile because no cutting is needed, which excludes bleeding and reduces the infection risk. Additionally, cryotherapy has no side effects, is very safe and leaves no scars.

Two versions

The CryoPen by Bentlon is available in 2 versions:

- CryoPen by Bentlon M
- CryoPen by Bentlon X

The CryoPen by Bentlon M has 2 micro-applicators for treatments of 2-5 mm and 4-8 mm, targeting the most frequent size of warts. With the 8 grams gas cartridge filled with nitrogen oxide it is possible to treat more clients.

The CryoPen by Bentlon X offers optimal flexibility because of the 4 different micro-applicators to choose from, varying from 1-3 mm, 2-5 mm, 4-8 mm to 7-15 mm. this enables you to choose the exact right applicator for each wart. Also, you will have the choice between 2 gas cartridges of either 8 or 16 grams nitrogen oxide. Depending on the number of clients you can choose the most economical cartridge to work with. Apart from standard cleaning, the CryoPen by Bentlon does not require extra maintainance.

Podo equipment	CryoPen by Bentlon M	CryoPen by Bentlon X
Ref.no.	65300	65310
Vloeistof	Nitic oxide	Nitic oxide
Freezing power	1 mm per 5 sec.	1 mm per 5 sec.
Temperature	-89°C	-89°C
Cartridge	8 gram (100-120 sec.)	8 gram (100-120 sec.) 16 gram (220-240 sec.)
Applicators	White: 2-5 mm	Blue: 1-3 mm
	Green: 4-8 mm	White: 2-5 mm
		Green: 4-8 mm Yellow: 7-15 mm
Accessoires	Cartridge 8 gram (24 stuks)	Cartridge 8 gram (24 stuks)
	Applicators	Cartridge 16 gram (6 stuks) Applicators





Podiatry accessories



1. Bentlon Cabinet

The Cabinet is a highly stable mobile pedicure cabinet for equipment and accessories. It has a separate shelf on the side for (pedicure) equipment and the top section offers plenty of space for all the accessories and equipment you need. It has three very spacious drawers. The Cabinet is also available without a separate shelf on the side.

2. Foot control with cable

There is a connection for this easy-to-use foot switch on all Bentlon pedicure units, except the Gold series. You can use it to operate the start and stop function.

3. Wireless Foot Control Gold

You can order this wireless foot control with the three Bentlon Podo S Gold pedicure units. You can switch the unit, the motor and the spray on and off with it and reverse the direction of the motor.

4. Omnicut hard skin file

This attachment is only available for the Pododent 2 and 2L hand attachments and is used to file away the hard skin. It saves about 15 minutes per client as compared with cutting hard skin away with a knife.

5. Contra-angle

This attachment with built-in light source is only available for the Pododent 2 and 2L hand attachments. You can increase its speed up to 200,000 revs/min by a right-angled transfer.

6. Pana Spray

Pedicure motors with removable attachments must be lubricated after every 12-15 treatments. Pana Spray is the tried-and-tested means for this worldwide. A connecting piece is supplied with Bentlon's detachable motors so that the spray nozzle of the spray can fits perfectly.

7. Leg cushion

A soft cushion for a comfortable support of the lower leg.

8. Aluminium tray

Aluminium tray for instruments. Fits on the top of the Bentlon Podo S Platinum unit.

9. Aluminium instrument shelf

Shelf for the aluminium tray to keep the instruments in place.

10. Plastic disposable tray

Instrument tray with pre-shaped rack and division that can be disposed of after use.

11. Massage head support

A special massage head support is ideal when giving massages while the client is lying on the stomach. Because of the round shape and hole in the middle, the client is comfortable in the lying position and can breathe without any obstruction. Available in the same colours as the chair cover.

Chair covers


Towelling chair covers for Bentlon Silver and Platinum models.

Protective sets


Loose set of covers for head support, arm supports and leg support on the Bentlon Silver and Platinum models.



1. Ref. 65625

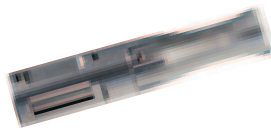

4. Ref. 66060

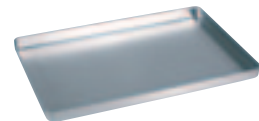

8. Ref. 66071



5. Ref. 66075




9. Ref. 66072



2. Ref. 60856



6. Ref. 66070

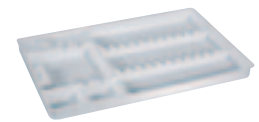

10. Ref. 66073



3. Ref. 66065



7. Ref. 66067



11. Ref. 61083.100 t/m 137

Website

Customerportal

All Bentlon customers can access the Dancohr dealer portal, where you will find everything you need to know about Dancohr products. View your order, look over contracts, update your customer loyalty card or even post requests/questions and request manuals for our products through the dealer portal.

Shared documents

But the dealer portal has much more to offer. You also have the possibility to download a variety of documents, from photographs and texts to logos, all of which provide you with active support in selling our products in your salon. You can even use these documents to design your own flyers, special offer campaigns and more.

News

Finally, the dealer portal also contains a news section, where you will find information on the latest developments regarding the various Dancohr brands, including Bentlon. Registration is not automatic. Go to www.dancohr.com to register for the dealer portal.

www.bentlon.com

In addition to the Bentlon image brochure, a Bentlon website has also been developed. The website contains descriptions of every product by category, accompanied by beautiful ambiance photographs.

Specifications & prices

The advantage of the Bentlon website is that it also contains specifications and prices for every product, convenient when looking for additional details on products or if you do not have the price list on hand.

News

When new products are launched on the market, they are always highlighted in the news section of the website. This section contains not only product introductions, but also other bits of news on Bentlon and Dancohr.

References

If you would like to experience Bentlon for yourself, you are welcome to visit the Dancohr Corporation in Weert or the Dancohr stand at the trade show! You can also look under the 'References' section on the Bentlon website, where you will find one or more Bentlon Specials for every province. This provides you with an overall picture of Bentlon in action!



Color swatches upholstery

	.100 White*
	.101 Sand
	.102 Pistachio
	.103 Mint
	.104 Aqua
	.105 Lemon
	.106 Yellow
	.107 Saffron
	.108 Orange
	.109 Topaz
	.110 Agave
	.111 Smaragd
	.112 Mais
	.113 Apricot
	.114 Mango
	.115 Terra
	.116 Fire
	.117 Cherry
	.118 Pearl
	.119 Silver
	.120 Chrome*
	.121 Stone
	.122 Dolphin
	.123 Laguna
	.124 Jeans*
	.125 Dark blue
	.126 Bordeaux
	.127 Ocean blue
	.129 Bright pink
	.130 Bright purple
	.131 Dark purple
	.132 Royal blue
	.133 Chocolate
	.134 Dark grey
	.135 Black
	.136 Beige
	.137 Creme

* .100/.120/.124 are standard. Other colours are on fee basis (see price list)



Want to know more about a specific treatment? Or how to clean a device properly? Or would you like to post a short video of a treatment on your website? YouTube has it all! The YouTube website contains instruction videos on various Bentlon treatments, techniques and cleaning methods. You can view them online or download them free.



This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

[illegible]



Dancohr Corporation BV
Marconilaan 8
6003 DD Weert
The Netherlands

Telephone General (+31) (0)495 434 600
Telephone Sales NL (+31) (0)495 434 601
Telephone Export (+31) (0)495 434 612

E-mail export@dancohr.com
Website www.bentlon.com

www.bentlon.com

Dealer/Distributeur

